


AndraTec *Advanced Interventional Technologies*

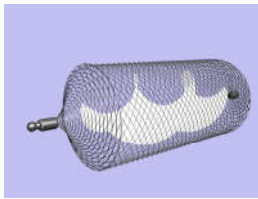

Pillow
AAA-Occluding Device

- Unique Endovascular peripheral occlusion device
- Highly Flexible and conformable within the vessel
- Low profile delivery system
- Optimal filling to promote fast thrombosis and prevent flow in variable locations
- The unique flexible and atraumatic design provides a barrier to blood flow thus enhancing vessel thrombosis
- Ease of use to ensure safe and atraumatic delivery and placement at the desired location
- Very low outward forces to prevent further pressure to the Aneurysmatic changes of the vessel

Technical Information:

Material:	Braided Nitinol
Diameter:	2.5 cm
Length:	8 cm
Profile:	8F sheath compatibility
Recommended Sheath:	8F / 90 cm
Recommended Guide Wire	0.035" (0.089 mm)
Recommended Sheath Introducer:	Ranging 8F

Pillow Occlusion Device Order Codes

Length Diameter	40mm	80mm
2.5 cm	PI-25-40	PI-25-80

SDS Pillow Delivery System Order Code

Length 160cm / 8F	SDS-1600
----------------------	----------

CE pending

Rev. I / 2019



AndraTec GmbH
Simmernerstr. 70
D-56075 Koblenz
Tel: +49 261 29675 172
Fax: +49 261 29675 173
www.andratec.com